



Reliance®..... N-2
Adhesion Booster

Reliance..... N-3
Protective Sealant

Reliance..... N-4
Bracket Adhesives Light Cure

Reliance..... N-7
Bracket Adhesives Chemical Cure

Reliance..... N-8
Cements

Reliance..... N-10
Specialty Procedures

Reliance..... N-11
Lingual Retention

Reliance..... N-12
Tooth Preparation

Bonding Auxiliaries N-13


SECTION - N

BONDING

Supplies



RELIANCE ADHESION BOOSTER

RELIANCE ADHESION BOOSTER

ITEM	DESCRIPTION	ITEM NUMBERS
	<p>Assure Plus All Surface Bonding Resin All Surface Light Cure Bonding Primer and Enhancer</p> <p>The Most Versatile Resin Ever Developed</p> <ul style="list-style-type: none">• Enhances bond strength on all surfaces• Tolerates moisture contamination• Primes and conditions any tooth surface for bonding• Effective on any enamel type• Effective on any artificial substrate• Eliminates hydrofluoric acid with porcelain crowns <p>Clinical Applications</p> <ul style="list-style-type: none">• Normal or atypical enamel• Dentin or cementum• Eliminates additional conditioners with:<ul style="list-style-type: none">• Stainless - Composite• Amalgam - Acrylic• Gold - Pontics• Eliminates hydrofluoric acid application with:<ul style="list-style-type: none">• Porcelain - Zirconia• Lithium Disilicate <p>1× bottle</p> <p>6 mL bottle</p> <p>2 mL bottle</p>	<p>PLUS</p> <p>PLUS2ML</p>
	<p>Assure Plus Unidose 50× 2 mL</p>	<p>UAP50</p>

RELIANCE PROTECTIVE SEALANT

RELIANCE PROTECTIVE SEALANT

ITEM	DESCRIPTION	ITEM NUMBERS
 <p data-bbox="183 777 375 818">Original Pro Seal cures in 360 to 420 range</p>	<p data-bbox="518 476 774 518">PRO SEAL® Light Cure Highly Filled Sealant</p> <p data-bbox="518 534 837 555">Enamel Protection From Start to Finish</p> <ul data-bbox="518 559 949 673" style="list-style-type: none"> • Two year surface protection • University tested in over 15 studies and clinically proven • Sets without oxygen inhibited layer - smooth, hard finish • Fluoresces under UV light • Patented catalyst system <p data-bbox="518 694 694 714">Clinical Applications</p> <ul data-bbox="518 721 981 808" style="list-style-type: none"> • Routine enamel protection • Hygiene challenged patients • Mitigates decalcification • Can be applied under bracket as a primer or around bracket <p data-bbox="518 839 614 880">PRO SEAL 6 mL bottle</p>	<p data-bbox="1356 859 1396 880">PRO</p>
 <p data-bbox="175 1274 391 1315">L.E.D. Pro Seal cures with any curing light, including L.E.D.</p>	<p data-bbox="518 963 774 1004">L.E.D. PRO SEAL® Light Cure Highly Filled Sealant</p> <p data-bbox="518 1036 614 1056">6 mL bottle</p>	<p data-bbox="1340 1025 1412 1046">LEDPRO</p>

RELIANCE BRACKET ADHESIVES - LIGHT CURE

RELIANCE BRACKET ADHESIVES - LIGHT CURE

ITEM	DESCRIPTION	ITEM NUMBERS
------	-------------	--------------



GoTo®

Light Cure Bracket Adhesive

GoTo® is a light cure adhesive with superior handling properties. It does not allow bracket flotation once placed on the tooth or during flash removal. A fraction of the retail cost of leading, competitive light cure pastes!

Newly Developed - Superior Working Properties


- Proven bond strength and consistency
- Unique viscosity easily penetrates bracket pads
- No bracket drift - stays in place until cured
- New proprietary microfiller construction
- One piece flash removal - "ropes" off bracket periphery
- Translucent finish matches any tooth shade
- Maximum bond strength
- Easy removal and cleanup
- UV disclosing agent - fluoresces on tooth

Clinical Applications

- Bracket and tube bonding
- Fixed lingual retainer bonding
- Bonds all metal, ceramic and composite attachments
- Now fluoresces under black light for complete cleanup!



Translucent

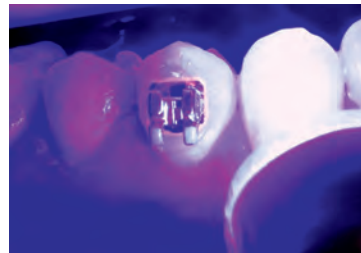

Blacklight Fluorescent

GoTo Light Cure Bracket Adhesive Push Syringe (1x 4 gm syringe)

GTP

GoTo Light Cure Bracket Adhesive Paste in Tips (10x 4 gm per tip)




GTT



For recommended primer Assure Plus, see page N-2.

RELIANCE BRACKET ADHESIVES - LIGHT CURE


RELIANCE BRACKET ADHESIVES - LIGHT CURE

ITEM	DESCRIPTION	ITEM NUMBERS
	<p>Light Bond™ Light Cure Bracket Adhesive</p> <p>Time Tested Performance</p> <ul style="list-style-type: none"> • Unrivaled bond strength • Highly filled paste • Smooth tacky viscosity • Consistent successful bonds • Easy flash and debonding cleanup • Designed for Doctor or staff bracket placement without drift <p>Clinical Applications</p> <ul style="list-style-type: none"> • Routine bracket and tube bonding • Intraoral bonding of metal, ceramic or composite attachments • Indirect bonding to stone/resin models <p>Light Bond Paste In Syringe (1× 5 gm syringe with fluoride)</p> <p>Regular Viscosity Syringe</p> <p>Medium Viscosity Syringe</p> <p>Thin Viscosity Syringe</p> <p>Blue Tinted Syringe</p>	
	<p>Light Bond Paste In Tips (11× 5 gm in tips with fluoride)</p> <p>Regular Viscosity Syringe</p> <p>Medium Viscosity Syringe</p> <p>Thin Viscosity Syringe</p> <p>Blue Tinted Syringe</p>	<p>LBPPF</p> <p>LBMPPF</p> <p>LBTPPF</p> <p>LBPPPF</p>
	<p>Light Bond Regular Kit (20 gm paste with fluoride, 2× 3cc resin, 23 gm liquid etchant, mixing pads, spatulas & brushes)</p> <p>With Sealant Resin</p> <p>With Assure Plus</p> <p>Light Bond Medium Kit (20 gm paste with fluoride, 2× 3cc resin, 23 gm liquid etchant, mixing pads, spatulas & brushes)</p>	<p>LBPTF</p> <p>LBMPTF</p> <p>LBTPTF</p> <p>LBPTF</p>
	<p>Push Syringe</p> <p>Push Syringe With Assure Bonding Resin</p>	<p>LBEMPF</p> <p>LBEMPF W-AS</p>

For recommended primer Assure Plus, see page N-2.

RELIANCE BRACKET ADHESIVES - LIGHT CURE

RELIANCE BRACKET ADHESIVES - LIGHT CURE

ITEM	DESCRIPTION	ITEM NUMBERS
	<p>Flow Tain™ Light Cure Flowable Microfil Composite Flow Tain™ Low Viscosity is a flowable light cure composite. The composite is highly polishable which will reduce wear and enhance patient comfort. When cured, Flow Tain Low Viscosity has low modulus of elasticity and is not as brittle as conventional composites. Even though Flow Tain is a flowable material, the thixotropic properties prevent slumping for better control. Using the tips provided, Flow Tain is easily injected onto lingual surfaces and across fiber splints. Comes in two different types of viscosities: Regular Viscosity and Low Viscosity.</p> <p>Clinical Applications</p> <ul style="list-style-type: none"> • Lingual Retainer Adhesive • Attachments on thermoplastic tooth aligning products • Splinting material (e.g fibers) • Indirect bonding adhesive • Archwire Stops • Clear aligner attachments <p>Regular Viscosity Paste</p> <p>Low Viscosity Paste 1.5 gm syringe with 10×19 gauge tips</p>	<p>FT FTLV</p>
	<p>19 Gauge Tips 10× tips</p>	<p>19TIPS10</p>

RELIANCE BRACKET ADHESIVES - CHEMICAL CURE

RELIANCE BRACKET ADHESIVES - CHEMICAL CURE

ITEM	DESCRIPTION	ITEM NUMBERS
	<p>Rely•A•Bond® No Mix Bracket Adhesive</p> <p>1 Step/Contact Adhesive</p> <ul style="list-style-type: none"> • Contact cure - primer to paste • Chemical cure - no light required • Saves time and material: <ul style="list-style-type: none"> » No mixing • Minimizes adhesive per bond <ul style="list-style-type: none"> • Highly filled paste: <ul style="list-style-type: none"> » Increased bond strength » Tacky viscosity » Prevents bracket drift • Choice of packaging: <ul style="list-style-type: none"> » Push syringe » Syringe tips 	
	<p>Rely•A•Bond Paste In Syringe Kit (14 gm paste, 16 gm primer, 9 gm etchant, spatulas, brushes and mixing pads)</p>	<p>RPF</p>
	<p>Rely•A•Bond Paste in Kit (14 gm paste, 16 gm primer, 9 gm etchant, spatulas, brushes and mixing pads)</p> <p>Separate Refill Components Rely•A•Bond Paste In Syringe (1× 3.5 gm push syringe)</p> <p>Rely•A•Bond Thin Paste In Syringe (1× 3.5 gm push syringe)</p> <p>Rely•A•Bond Paste In Tips (14 gm syringe tips)</p> <p>Rely•A•Bond Primer with Fluoride (1× 16 gm Bottle)</p>	<p>RTF</p> <p>1RPP</p> <p>1RPPT</p> <p>RPT</p> <p>RPRF</p>


RELIANCE CEMENTS

RELIANCE CEMENTS

ITEM	DESCRIPTION	ITEM NUMBERS
	<p>Ultra Band-Lok™ Light Cure Compomer Band Cement Industry Leader in Band Cement</p> <ul style="list-style-type: none"> • Single Paste System • Light Activated • Resin Base - Glass Ionomer • High Fluoride Release • Natural or Blue Shade <p>Bonds To Metal Ultra Band-Lok® will bond chemically to stainless steel bands. Removal is all in one piece with cement in the appliance, making cleanup effortless.</p> <p>All Purpose Band Cement Ultra Band-Lok® can also be used to bond large acrylic appliances and as an occlusal build-up material.</p>	
	<p>Clinical Applications</p> <ul style="list-style-type: none"> • Band cementation • Occlusal build ups • Acrylic appliances - Hyrax, MARA, RPE, etc. • Lab block out material • Non-fluoride upon request 	 <p>Simplified Cleanup</p>  <p>Occlusal Buildup</p>
	<p>Ultra Band-Lok Syringe Push Syringe 1× 5 gm syringe</p> <p>Push Syringe Blue Shade with Fluoride 1× 5 gm syringe</p> <p>Push Syringe Blue Shade without Fluoride 1× 5 gm syringe</p> <p>Ultra Band-Lok Luer Loc Syringe With Dispensing Tips 1× 5 gm syringe</p> <p>Blue Shade With Dispensing Tips 14 gm tips</p> <p>(5) Luer Loc Dispensing Tips Luer Loc Dispensing Tips 14 gm tips</p> <p>(50) Luer Loc Dispensing Tips 14 gm tips</p>	<p>UBLPP</p> <p>UBLBPP</p> <p>NFUBLBPP</p> <p>LLUBLP</p> <p>LLUBLBP</p> <p>TIPS50</p> <p>TIPS5</p>

RELIANCE CEMENTS

RELIANCE CEMENTS

ITEM	DESCRIPTION	ITEM NUMBERS
 <p>A)</p>	<p>GlassLok™ Dual Cure Glass Ionomer Band Cement Hydrophilic Powder/Liquid Cement</p> <ul style="list-style-type: none"> • Secure band adhesion • High compressive and tensile strength • Bonds in wet or dry fields • Resists washout from under bands • Contains fluoride – protects enamel • Bonds chemically to metal, reducing removal cleanup time • Dual cure - chemical or light activated • Natural or blue shade <p>GlassLok Powder A) Powder 15 gm</p> <p>Powder Economy Size 100 gm (1x powder)</p>	<p>GLP15 GLP100</p>
 <p>B)</p>	<p>GlassLok Liquid B) Liquid Natural Shade 15 gm</p> <p>Liquid Natural Shade Economy Size 25 gm</p> <p>Liquid Blue Shade 15 gm</p> <p>Liquid Blue Shade 25 gm (1x bottle)</p>	<p>GLL15 GLL25 GLL15B GLL25B</p>
 <p>C)</p>	<p>GlassLok Starter Kit Natural Shade</p> <p>Blue Shade (1x 15 gm powder, spatulas, measuring scoop, 1x 15 cc liquid, mixing pads)</p> <p>GlassLok Economy Kit C) Natural Shade</p> <p>Blue Shade (1x 100 gm powder, spatulas, measuring scoop, 3x 25 gm liquid, mixing pads)</p>	<p>GLS GLSB GLE GLEB</p>

RELIANCE SPECIALTY PROCEDURES

RELIANCE SPECIALTY PROCEDURES

ITEM	DESCRIPTION	ITEM NUMBERS
------	-------------	--------------



Bond Aligner®
Light Cure Aligner Adhesive
Bonds Attachments to Clear Aligners

- Secure aesthetic bonds
- Bond to any clear aligner material
- No primers or conditioners
- Clear light cure paste
- High modulus of elasticity:
 - Conforms to aligner distortions
 - Secure bonds throughout
- Bond ceramic, metal or plastic attachments

Bond Aligner Syringe (1x 1.4 gm syringe, 5x tips)

BA



C&C™ Composite and Chain System
Molded Base and Chain Extrusion System
Unparalleled Strength

The C&C System does not use protruding metal attachments. Our patent pending design molds a base with composite - adapting flush to any palatal or buccal exposure.

Painless Removal

Instead of painfully breaking a bracket off a mobile tooth - the C&C™ is ground down with a bur - resulting in less patient discomfort.

Disposable Attachment Mold with Chain

Large (1x pad)


DCCP

Small (1x pad)

DCCO


1.5 g FlowTain™ Flowable Composite

FT

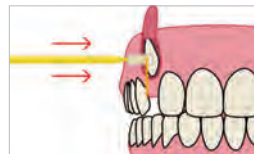

RELIANCE LINGUAL RETENTION

RELIANCE LINGUAL RETENTION

ITEM	DESCRIPTION	ITEM NUMBERS
	<p>Reliance Ortho FlexTech® Flexible Lingual Retainer - A New Generation in Fixed Retention</p> <ul style="list-style-type: none"> • Secure and stable long term passive retention • Patented interlocking chain design • Adapts to lingual anatomy – No bending or burnishing • No models or impressions • Flexes interproximally – Resists breakage and bond failures • No reciprocal forces – allows for minor biological movement – will not take an active bend • Stainless steel offers increased tensile and torsional strength • Low profile – patient comfort <p>Stainless Steel Spool (30" spool)</p> <p>White Gold Spool (30" spool)</p> <p>Nickel Free Ortho FlexTech (Retention System 30" Length)</p>	<p>SSOFT</p> <p>OFT2</p> <p>NFOFT2</p>
	<p>Reliance RETAINIUM® XL Optimum Strength Beta Titanium Retainer Wire</p> <ul style="list-style-type: none"> • Fixed Retention with Titanium Molybdenum • Flat .027" x .016" ribbon arch • Low profile – comfortable and functional • Increased strength • MRI compatible • Easy adaptation • Nickel free titanium construction <p>Retainium XI Titanium Retainer Wire (10x 6" wires)</p>	<p>XLRET</p>
	<p>Reliance RETAINIUM Original TM Alloy Retainer Wire</p> <ul style="list-style-type: none"> • Lingual Titanium Molybdenum Retention • Lowest profile • .027" x .011" ribbon arch • Adaptability • Flexes Interproximally • Nickel free • MRI compatible <p>Retainium Alloy Retainer Wire (10x 6" wires)</p>	<p>RET</p>
	<p>Reliance® Bond•A•Braid® Lingual Retainer Wire</p> <p>The dead soft wire easily adapts and is excellent for semi-permanent orthodontic splinting of the upper and lower incisors, diastema maintenance and other tooth splinting.</p> <ul style="list-style-type: none"> • Lingual retainer • .027" x .011" ribbon arch wire • 8-strand braided wire • Dead soft – adapts easily <p>Bond•A•Braid Lingual Retainer Wire (10x 6" wires)</p>	<p>BAB2</p>



RELIANCE TOOTH PREPARATION

RELIANCE TOOTH PREPARATION



ITEM	DESCRIPTION	ITEM NUMBERS
	<p>Blue Gel Etchant</p> <ul style="list-style-type: none"> Controlled and even application 37% phosphoric acid Excellent for lingual retention or large acrylic appliances Contains antimicrobial B. A. C. <p>Gel Etchant (1× 5 gm syringe, 10× tips)</p>	<p>BEGS</p>
	<p>Gel Etchant</p> <ul style="list-style-type: none"> 37% phosphoric acid Controlled application <p>Gel Etchant Syringe (1× syringe, 20× tips)</p> <p>Gel Etchant Bottle (1× 9 gm bottle)</p> <p>Gel Etchant JAR (1× 18 gm jar)</p>	<p>EGSS</p> <p>EG9</p> <p>EG18</p>
	<p>Liquid Etchant</p> <ul style="list-style-type: none"> 37% phosphoric acid Green disclosing tint Regular or cherry flavor <p>Liquid Etchant Small (1× 9 gm bottle)</p> <p>Liquid Etchant Large (1× 23 gm bottle)</p> <p>Liquid Etchant Large Cherry Flavor (1× 23 gm bottle)</p>	<p>EL9</p> <p>EL24</p> <p>EL24C</p>
	<p>Porc Etch™ Porcelain Etching Agent with Barrier Gel</p> <ul style="list-style-type: none"> Preps porcelain crowns for bonding <p>Porc Etch with Barrier Gel (1× 7 gm bottle, 1× 3 mL bottle)</p>	<p>POE</p>
	<p>Porcelain Conditioner For Bonding to Ceramic Crowns</p> <ul style="list-style-type: none"> Effective with: porcelain, zirconia and lithium disilicate <p>Porcelain Conditioner (1× 8 cc bottle)</p>	<p>PO</p>

RELIANCE TOOTH PREP / BONDING AUXILIARIES

RELIANCE TOOTH PREPARATION






ITEM	DESCRIPTION	ITEM NUMBERS
	<p>Assure The Original Universal Bonding Resin</p> <ul style="list-style-type: none"> • Bonding primer and enhancer • Typical or atypical enamel • Metal and composite restorations <p>Assure Universal Bonding Resin Large (1x 6 mL bottle)</p>	AS
	<p>Light Bond Sealant Light Cure Resin</p> <ul style="list-style-type: none"> • Bonding primer and sealant • Available with or without quartz filler <p>Regular Sealent (1x 3 gm bottle)</p> <p>Filled Sealent (1x 3 gm bottle)</p> <p>Regular Sealent (1x 7 gm bottle)</p> <p>Filled Sealent Item (1x 7 gm bottle)</p>	<p>LBS3F</p> <p>LBF5F</p> <p>LBS7F</p> <p>LBF7F</p>

BONDING AUXILIARIES

ITEM	DESCRIPTION	ITEM NUMBERS
	<p>Mini Sponge Applicators (200 Sponges in Plastic Jar)</p> <p>Specially designed for application of liquids, etches, and resins to teeth. Using tweezers, each sponge can be picked out from premoistened jar and dipped in liquid before applying to desired area. Mini sponges may be used during bonding to apply etchant or any liquid, such as primers, conditioners, etc.</p>	480-015
	<p>Cheek Retractors</p> <p>For treatment with teeth in open or closed bite. Self-retaining. Can be cold sterilized. Excellent for photography.</p> <p>Adult Cheek Retractors (2 per pack)</p> <p>Teen Cheek Retractors (2 per pack)</p> <p>Pedo Cheek Retractors (2 per pack)</p> <p>Assortment Pack (one of each size)</p>	<p>300-020</p> <p>300-022</p> <p>300-021</p> <p>300-023</p>

BONDING AUXILIARIES

BONDING AUXILIARIES

ITEM	DESCRIPTION	ITEM NUMBERS
	<p>Lip Retractors Transparent plastic lip retractors have long handles so they can be held easily. Excellent for photography. Rounded edges. Can be cold sterilized. (2 per pack)</p> <p>Standard Lip Retractors</p> <p>Mini Lip Retractors</p>	<p>300-025</p> <p>300-026</p>
	<p>Retract EEZ Photo Lip Retractor Excellent for occlusal photographs with easy to grip handles, a great anatomical fit, plus angled handles. Can be cold sterilized, as well as autoclaved, though heat sterilization may cause hazing of the material.</p>	<p>300-029</p>
	<p>Tongue Shield – Latex Free Combination of two bite blocks and a vertical shield. The shield prevents the tongue from touching the lingual surface of the anterior teeth. Center hole allows for insertion of saliva ejector. Excellent during bonding procedures. Can be cold sterilized. (Pink, 3 per pack)</p> <p>Adult/Teen</p> <p>Child/Pedo</p>	<p>300-030</p> <p>300-031</p>
	<p>Nola™ Dry Field System The Nola Dry Field system completely eliminates saliva and fully exposes both dental arches making it ideal for all bonding applications. The system is easy to assemble and can be inserted and used by a single operator. Autoclavable.</p> <p>Adult</p> <p>Pedo</p> <p>Connector</p>	<p>330-029</p> <p>330-030</p> <p>330-031</p>
	<p>Deluxe Adjustable Cheek Retractors Unique bow design allows adjustment to ensure comfortable fit. Can withstand temperatures up to 400°F. Fully sterilizable. (1 per pack)</p>	<p>300-015</p>